



SKP011

SECTIONAL MEMORY
FOAM MATTRESS



SAIKANG / Medical

Healthier
your life

SKP011

FEATURES AND HIGHLIGHTS

Technical parameters

▪ Size (LxWxH) 1925×900×140mm

Zippered Cover

▪ Material PU
▪ Color grey or dark blue
▪ Size (LxWxH) 1925×900×140mm

Upper Layer

▪ Material visco-elastic memory foam
▪ Color blue, green, black
▪ Thickness 40mm
▪ Density 42-45kg/m³

Firm Base Layer

▪ Material polyurethane foam
▪ Color creamy-white
▪ Thickness 100mm
▪ Density 40kg/m³


Upper Layer

a) High density: 42-45kg/m³
b) Unique thermo active properties to provide good support and excellent pressure relief



Firm Base Layer

a) High density: 40kg/m³
b) Very elastic and thermosetting. Sectional design works to make the mattress and bed surface fit well in different positions such as back-rest and knee-rest adjustment.

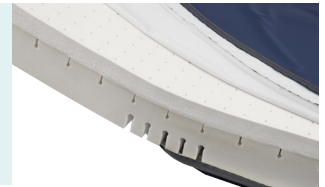



Zippered Cover

Air through, water proof, soft, wind-resistant, strong tensile strength and anti-infection. Long zipper design, easy to remove and wash.



Size

Customized size acceptable

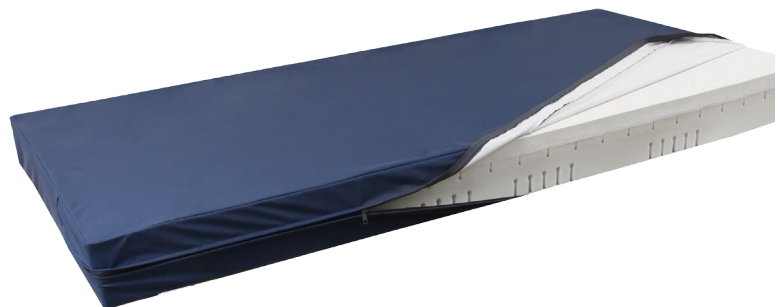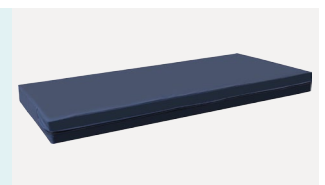



CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue