



B8e

MULTIFUNCTIONAL ICU ELECTRIC BED

Registration code: SKD-C




SAIKANG / Medical

Healthier
your life


B8e

FEATURES AND HIGHLIGHTS

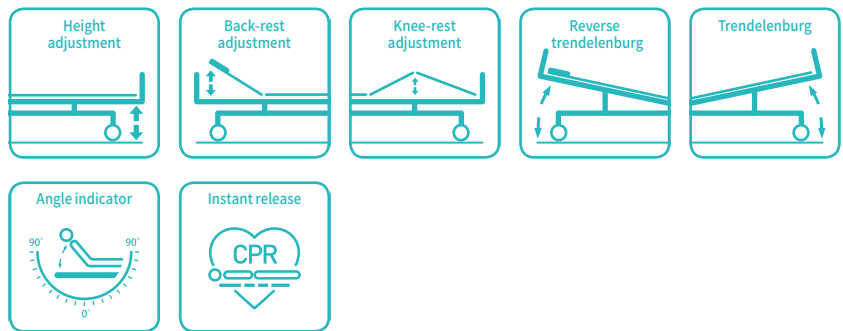

Technical parameters

▪ External size (LxW):	2170x1065 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	435-725 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse Trendelenburg:	0-12°
▪ Trendelenburg:	0-12°
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	

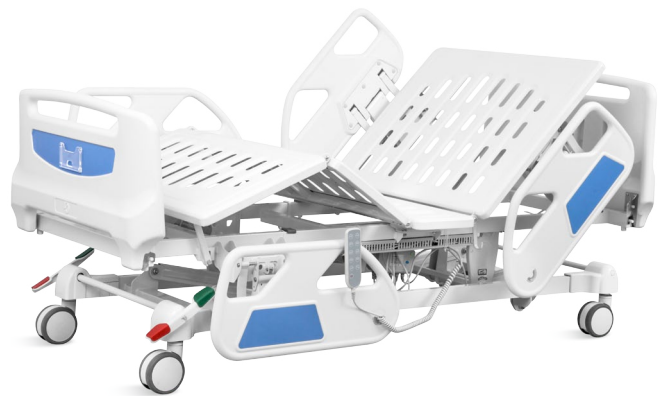


Technical configuration

<input checked="" type="checkbox"/> Linear Motor	4 pcs
<input checked="" type="checkbox"/> 5" Double Side Castors, Central Locking	4 pcs
<input checked="" type="checkbox"/> ABS Side Rail Extra Height	1 set
<input checked="" type="checkbox"/> ABS Bed Ends with Integrated Control Panel	1 set
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	2 pcs
<input checked="" type="checkbox"/> Central Locking Pedal	2 pcs
<input checked="" type="checkbox"/> Electric CPR	




Optional

<input type="checkbox"/> Manual CPR	
<input type="checkbox"/> Weighing System	
<input type="checkbox"/> Back-up Battery	
<input type="checkbox"/> X-ray Bed Surface	
<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan Traction Frame	
<input type="checkbox"/> Directional Wheel	

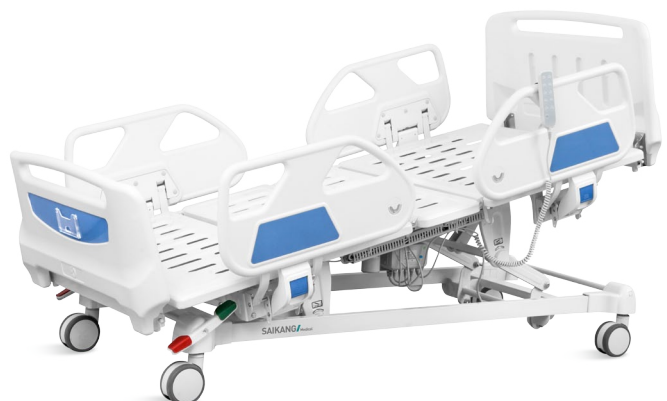

Design

E-shaped large guardrail, with angled display and built-in controller, meets safety standards, providing comforting feeling for patients.

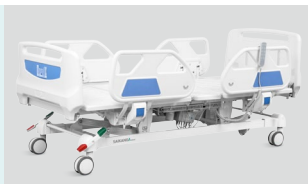**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth wedding. Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**


Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Angle Indicator**

Patented angle indicator.

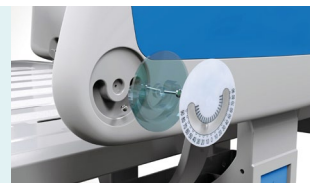**Wheels**

Secure® Castors
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Durable Brakes**

Break spring keep qualified after 10,000 times uninteruped break.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Side Rail**

Full ABS material, tuck away side rail with built-in controller. From top of the guardrail to bed surface is about 345mm, more height provides better protection for patient. Bottom can endure 100,000times press, buffer, built-in angle indicator.

**Head & Foot Board**

Full ABS material, Bottom can endure 100,000times press
Head and foot board can be reversable.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue