



# GT3k

## MANUAL ORTHOPEDIC TRACTION BED

Registration code: SKD-4C




**SAIKANG** Medical

Healthier  
your life

# GT3k

## FEATURES AND HIGHLIGHTS

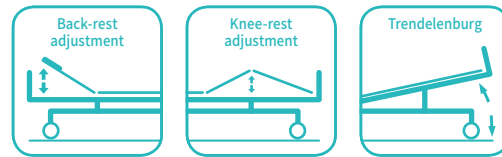

### Technical parameters

▪ External size (LxWxH)	2270x970x500 mm
▪ Mattress platform	1925x830 mm
▪ Safe working load	250 kg
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)
▪ Trendelenburg	0-6°



### Technical configuration

<input checked="" type="checkbox"/> Manual cranks	3 pcs
<input checked="" type="checkbox"/> 5" double side castors, central locking	4 pcs
<input checked="" type="checkbox"/> Central locking pedal	2 pcs
<input checked="" type="checkbox"/> Aluminium alloy ABS head and foot board with safe lock	1 set
<input checked="" type="checkbox"/> Aluminium alloy foldable side rail	1 set
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	2 pcs
<input checked="" type="checkbox"/> Aluminium balkan frame	1 set



**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth welding Solid bed base made of 5.0 mm stamped rotating plate.

**Balkan Frame**

Aluminum Balkan Frame, which can adjust the longitudinal, lateral and height positions in all directions, and can accurately adjust the most comfortable posture of the patient.

**Hand Control**

Retractable hand crank, 100,000 times remaining intact  
Extreme protection device, sounds when tipped to the max position,  
Strong pallets to increase the bed.

**Bumper Cover**

Reduce the impact of the collision, safe and slinky, protect the bed and patient during transportation.

**Platform**


Stamping forming platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform.

**Guardrail**

Aluminum alloy side rail, 6-staff, collapsed type with safe lock.

**Head & foot board**

Aluminium alloy flat tube, MDF material fire-proof plate. size of the head&foot board can be customized.

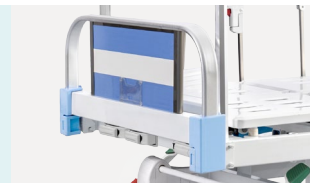**Wheels**

Castors with TPR tire, no worn-out after running 30KM, united forming without bolts.  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time.

**Brake**

Three position braking system, transporting can maintain a straight line.  
ABS cover on castors joint, prevent the dust or hair stuck in. Easy to clean.

**Anti-slip bed feet cover**

Non-skid, more stable in the ground.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue