



A99-5

# ELECTRIC OBSTETRIC BED

Registration code: SKF-C

**SAIKANG** / Medical

Healthier  
your life

# A99-5

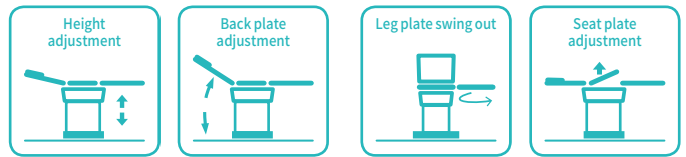
## FEATURES AND HIGHLIGHTS

### Technical parameters

▪ External (LxW):	1880x600mm
▪ Table height:	680-940mm
▪ Seat plate up folding:	0-20°
▪ Seat plate down folding:	0-9°
▪ Back plate up folding:	0-45°
▪ Back plate down folding:	0-10°
▪ Leg plate abduction:	0-90°
▪ Load bearing:	≤300kg

### Technical configuration

<input checked="" type="checkbox"/> Arm support	2 pcs
<input checked="" type="checkbox"/> Handle	2 pcs
<input checked="" type="checkbox"/> Leg support	2 pcs
<input checked="" type="checkbox"/> Stirrup	2 pcs
<input checked="" type="checkbox"/> Slider	4 pcs
<input checked="" type="checkbox"/> Tightener	2 pcs
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> Foot switch	1 pc
<input checked="" type="checkbox"/> Stainless steel basin	1 pc



**Material**

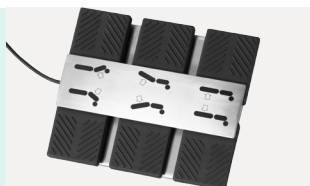
The chassis cover is made of high quality 304 stainless steel, which is welded by the robot seamlessly. The shape is beautiful and stable.

**Mattress**

Polyurethane foam mattress, waterproof and seamless, easy cleaning and disinfection.

**Foot pedal controller**

All the functions are controlled by hand controller and foot switch.

**Braking system**

The pedal controls the brake function of the casters.

**Control System**

Electric control the height lifting, seat plate up and down, back plate up and down, convenient operations.

**Arm Support**

The arm support could be adjusted into multiple angles with the straps, it can be used to comfortably fix the puerpera's arm at various angles.

**Leg supports**

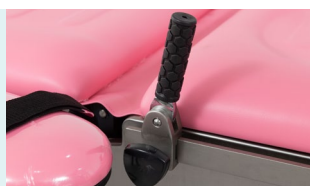
The leg support is integrated with stirrup function, and the delivery bed is equipped with a pair of stirrups.

**Auxiliary Table**

The auxiliary table can be folded 90° outwards, used as a instrument table, detachable, and adapt to a variety of surgical environments.

**Pull handles**

Pull handles can be folded to horizontal position. Easy patient mobilisation.

**Sewage basin**

Equipped with stainless steel sewage basin and other various accessories.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue