



Comfort series

# B4w

## MANUAL BED


Registration code: SK-CD1



**SAIKANG** / Medical

Healthier  
your life

## B4w Manual bed

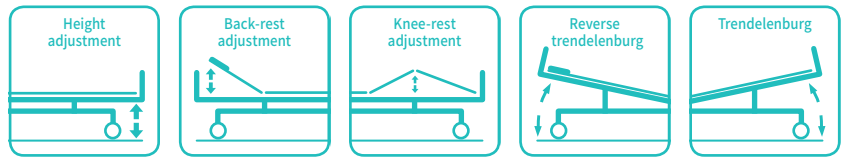

### Technical parameters

▪ Size (LxW):	2180x960 mm
▪ Bed height adjustment:	460-710 mm
▪ Mattress platform:	1925x830 mm
▪ Safe working load:	250 kg
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse trendelenburg:	0-12°
▪ Trendelenburg:	0-12°

### Technical configuration

<input checked="" type="checkbox"/> Manual crank	4 pcs.
<input checked="" type="checkbox"/> 5" hardshell castors with break	4 pcs.
<input checked="" type="checkbox"/> Foldable side rail with ProErgo-Handle™	1 set
<input checked="" type="checkbox"/> ABS bed ends with safe lock	1 set
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs.
<input checked="" type="checkbox"/> Drainage hook	2 pcs.

### Optional

<input type="checkbox"/> Directional wheels function	
<input type="checkbox"/> Monkey pole	
<input type="checkbox"/> Aluminium balkan frame	

### Design

Patented L-shaped foldable side guardrail: bended tube, solid installing, no shaking, anti-pinch, protect patient's safety.



### Hand control

Retractable hand crank, 100,000 times remaining intact. Extreme protection device, sounds when tipped to the max position.



### Bed surface

Stamped forming platform, smooth surface, no welding scar, safe working load 250 kg. Optional: ABS platform and resin X-ray translucent platform.



### Painting

Baosteel from Fortune Global 500. 11 process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, painting can resist 50 kg impact.



### TPR tire castors

TPR tire and save anti-winding hard shell, united forming without bolts. Pass Dynamic test: bearing 120 kg run 30 km and pass obstacles 500 times.



### Virgin ABS plastic

Environmental plastic can be sterilized in 100°C. Tensile upto 30MPa.



**SAIKANG** Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China  
+86 18021231901 · 24/7  
sales@chinasaikang.com  
www.saikangmedical.com



Company video



General catalogue