



# X9x

## MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C



**SAIKANG** / Medical


Healthier  
your life

# X9x

## FEATURES AND HIGHLIGHTS

### Technical parameters

▪ External size (LxW):	(2040-2240)x1020 mm
▪ Mattress platform:	(1820-2020)x820 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	500-950 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-20° (±10°)
▪ Reverse Trendelenburg:	0-18°
▪ Trendelenburg:	0-18°
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	

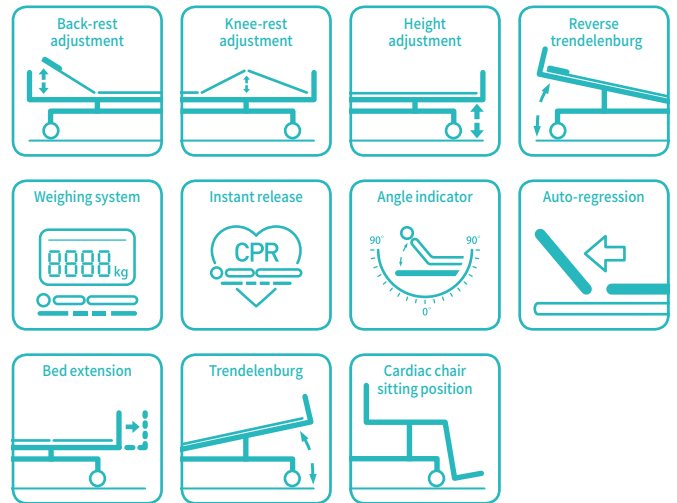

### Technical configuration

<input checked="" type="checkbox"/> Linear LINAK Motor	6 pcs
<input checked="" type="checkbox"/> 6" Double Side Castors, Central Locking	4 pcs
<input checked="" type="checkbox"/> ABS Side Rail with Integrated Control Panel	1 set
<input checked="" type="checkbox"/> ABS Bed Ends One Step Removable	1 set
<input checked="" type="checkbox"/> Ø75mm Wall Bumper Wheel	4 pcs
<input checked="" type="checkbox"/> IV Pole	1 set
<input checked="" type="checkbox"/> IV Pole Prevision	4 pcs
<input checked="" type="checkbox"/> Drainage Hook	4 pcs
<input checked="" type="checkbox"/> 5th Directional Wheel Pedal	2 pcs
<input checked="" type="checkbox"/> Central Locking Pedal	2 pcs
<input checked="" type="checkbox"/> X-ray Translucent Mattress Platform	1 set
<input checked="" type="checkbox"/> Back-up Battery	1 set
<input checked="" type="checkbox"/> Night Light	1 set
<input checked="" type="checkbox"/> Bed Push Handle for Emergency Transfer	1 set
<input checked="" type="checkbox"/> Manual CPR	
<input checked="" type="checkbox"/> Electric CPR	



### Optional


<input type="checkbox"/> Weighing System With LCD Screen
<input type="checkbox"/> Data Storage System
<input type="checkbox"/> Monkey Pole

**Design**

Extra Long and height head foot board, with curve surrounded design, provide conforting feeling for patients.

**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

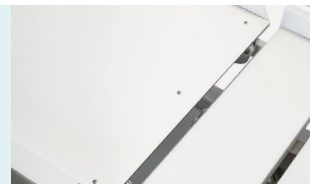**Welding | Safety**

Panasonic Robotic ensure, 360° full smooth wedding. IEC606011 Safty Standard.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**Electric motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**7 Inch touch screen**

Replacing the original button screen, more convenient operation. There are two language options, chinese and english. Screen can store the patient and doctor's information, users manuals, weight record

**Brake**

Break spring keep qualified after 10,000times uninteruped break.

**Wheels**

TPR tire no worn out after running 30 km, save anti-winding hard shell, united forming without bolts. Pass dynamic test: Bearing 120 kg, run 30 km and pass obstacles 500 time.

**Bed surface**

X-Ray Translucent Mattress Platform.

**Side rail**

Tuck away side rail with buffer Built-in controller and angle indicator.

**ABS bed ends**

ABS bed ends One Step Removable, meet the first-aid requirements of high-end wards and ICU room.

**Bumper wheels**

Bumper wheel ,  $\phi 74$ mm, max. protect the bed and patient during transportation.

**Battery**

Battery - 18times movement with 100kg bearing, 130times movement without bearing.

**IV pole**

Stainless Steel Column with plastic hook.

**Bottom**

The shelf at the bottom of the bed can hold patients items, like shoes and something others.

**Lifting pole**

Lifting pole - Bearing 765N force.





## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue