



Oasis Series

Y6y8c

Electric Bed
SKD-C





SAIKANG / Medical

LEADING BRAND OF CHINA MEDICAL FURNITURE
Stock Code: 870098

QUALITY CONTROL

GREAT MULTI-CARE

Respects easy patient mobilisation:
One-button for patient mobilisation position -back plate raised, double-hole design on back part guardrail, which supports patient self standing.

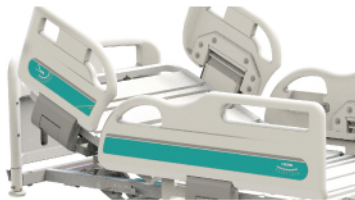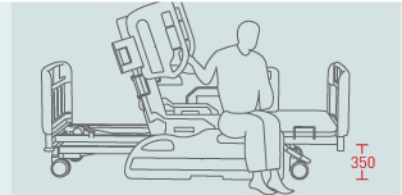


SAFETY GUARDRAIL

The designed guardrail fully considered to provide enough space and hand grasp for patient to get out of bed, and provide extremely high safety.
Two steps unlocking system, to prevent patients from accidentally lowering the guardrail (patented design).

DOUBLE REGRESSION FUNCTION

Auto-Regression 10+6 cm
When repositioning the backrest and thigh rest, the Y6Y8C makes the back plate retreat 10cm and the thigh plate retreats 6cm so to minimises the pressure and shear forces in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers.

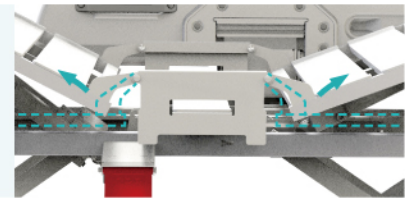


IEC60601

SAFETY

Follow IEC60601 Safety Standard
One-button emergency stop switch can quickly stop all going movement.

CPR

CPR emergency reset handle to immediately flatten the bed back part; located on both sides of the bed.
Electric CPR button is setted inside control panel, one button to make the back plate flat.

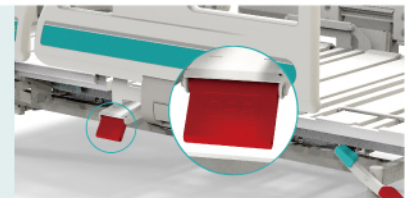

10,000_{times}

BRAKE

Three position braking system, transporting can maintain a straight line

ELECTRIC SYSTEM

American motor, with UL/ROHS/EN approval. Noise under 50dB, with lifetime 20,000times movement.



UL / ROHS / EN

50... | 20,000_{times}



BACK REST 30DEGREE POSITION


Help nurse to observe the bed position from distance, many research shows 30degree position will max. reduce the chance of lung infection.

RICH PARTS



Castors with TPR tire, no worn-out after running 30KM, united forming without bolts

Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 times.

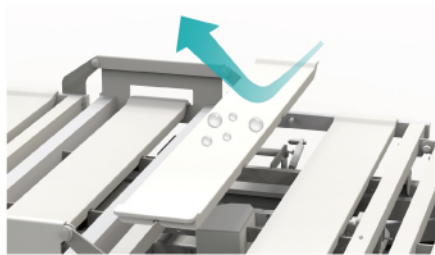

JOINT-FREE TECHNOLOGY


Removable plastic covers formed from 9 removable blow moulded lamella to minimise the slots where unwanted infections may hide.




Patent angle indicator, clearly show back-rest and Trendelenburg position



Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable

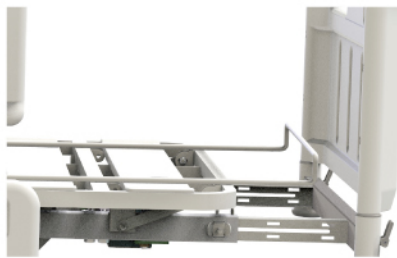

Extended bed platform up to 20CM, comfortable for the tallest of patients (option)

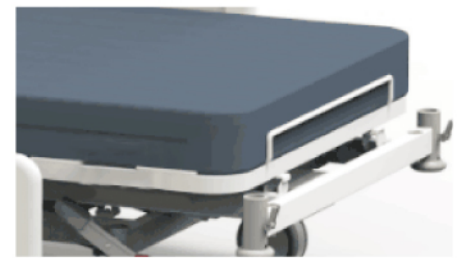

Mattress Holder.
Waterproof, antibacterial mattress.



Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall



5th and 6th extra side rail, for 100% fully protection



Socket on side rail for placing the patient cellphone or glasses



Plastic drainage hook, movable



Stainless Steel Column with plastic hook

DESIGN CONCEPT



Patient

Ergoframe® against pressure ulcers
Easy mobilisation
Available with a number of controls



Staff

Soft drop function
Easy clean concept
Safety features such as the Extender or Protector



Management


Low TCO – total cost of ownership
Universal bed for a low price
Overload protection



Back-rest



Knee-rest

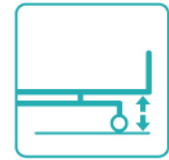


Height Adjustment



CPR



Bed Extension

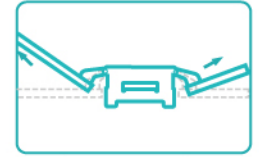

Auto-regression

TECHNICAL PARAMETERS

Bed Dimensions L*W*H	2130*1050*350-800mm
Bed Platform	1925*900mm
Back-rest Adjustment	0-75°(±10°)
Knee-rest Adjustment	0-35°(±10°)
Mattress Size	1925*880*10mm
Bumper Wheel	4*Φ80mm
Packing Size	2120*1080*450mm

TECHNICAL CONFIGURATION

American Motor	4pcs
5" Double Side <i>Secure</i> ® Castors	4pcs
ABS <i>Oasis</i> Side Rail	1set
ABS <i>Oasis</i> Bed Ends with safe lock	1set
Φ80mm Wall Bumper Wheel	4pcs
IV Pole	1set
IV Pole Prevision	4pcs
Drainage Hook	2pcs
Central Locking Pedal	1pc





INNOVATION OF SAIKANG

LEADING BRAND OF CHINA MEDICAL FURNITURE



Design Innovation



Research Innovation



Manufacturing Innovation



Maintenance Innovation



Service Innovation

JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD.

24hrs: 0086-18021231901 
www.saikangmedical.com

SAIKANG (China Headquarter)

No.35 Lehong Road,Zhangjiagang City,
Jiangsu Province,China.
TEL: +86-512-58525850
FAX: +86-512-58520713
E-mail: inquiry@jssaikang.com

SAIKANG (USA Branch)

Saikang Medical Technology, Inc
17800 Castleton St Ste 180,
City of Industry, Ca 91748
Tel: +1 949648-2776

SAIKANG (UK Branch)

Uk Healthward International Limited
Owner(s) Address DEPT. 111,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +44 2032399738

SAIKANG (East Africa Branch)

Saikang Medical Solutions (East Africa) Ltd.
Viraj Complex Unit 3 Mombasa Road,
P.O Box 14262-00800, Nairobi.
TEL: +254 733-704 690

SAIKANG (Mid East Branch)

Saikang Medical Supplies (Mid East) LLC
X22-S03 England Cluster,
International City, Dubai UAE
TEL: +20 128 8962497



Contact



Website

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.