



V2w5c

MANUAL BED

Registration code: SK-C1




SAIKANG Medical

Healthier
your life

V2w5c

FEATURES AND HIGHLIGHTS

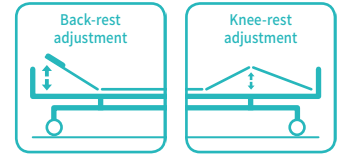

Technical parameters

▪ External size (LxWxH):	2130x1020x500 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	200 kg
▪ Back-rest adjustment:	0-70° (±5°)
▪ Knee-rest adjustment:	0-38° (±5°)

Technical configuration

<input checked="" type="checkbox"/> Manual Crank	2 pcs
<input checked="" type="checkbox"/> 5" Covered Castors	4 pcs
<input checked="" type="checkbox"/> L-shaped foldable side rail	1 set
<input checked="" type="checkbox"/> PP bed ends with safe lock with $\Phi 65$ bumper wheels	1 set
<input checked="" type="checkbox"/> PP bed platform with integrated mattress retainer	1 set
<input checked="" type="checkbox"/> IV pole 4 hooks	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs
<input checked="" type="checkbox"/> Drainage hook	2 pcs

Optional

- Auto-regression (back section)
- Trapeze bracket/pole

Optional color (MOQ100pcs- color changing)



QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



Design

The unique germ-free design meets hospitals' strict infection control requirement. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.



Steels | welding

Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



Plastic | Painting

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.
11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

Accessories



Platform

Integrated PP bed platform, adopting advanced manufacturing technology. 4-5 section platform, enables independent installation and easy maintenance.

Enlarge the seat area when back part is raised, providing more comforts for the patients.

Reduce the depression of the sponge at the seat panel.



Castors

5" covered casters. TPR tire no worn out after running 10 km.

Anti-winding hard shell cover, united forming without bolts.

Pass dynamic test: pass obstacles 40000 times, under the standard 200kg load of the whole bed.

Head & Foot Board

LG Chem fresh PP material, Plug in type. Equipped with the bumpers wheel, diameter 65mm, preventing collision.




Wall bumper wheel

Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall.



Guardrail

L-shaped foldable guardrail (Patent), bended tube, solid installing, 0 shaking, anti-pinch, ensure patient's safety.

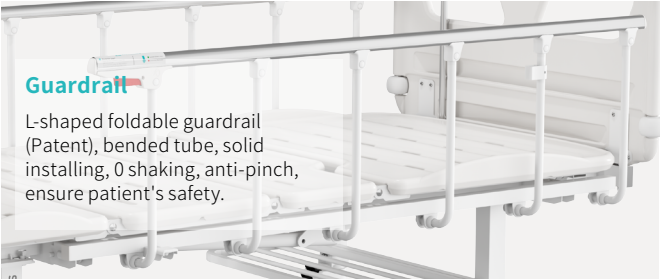

Mattress Retainer

Platform integrated mattress retainer at the both sides of the back and leg section, with a boarder of 33mm height.



IV Pole

SS IV pole, bearing 15kg.




Drainage Hooks

Mobile plastic drainage hooks.



Utility Shelf

Detachable utility shelf for patients' personal stuff.

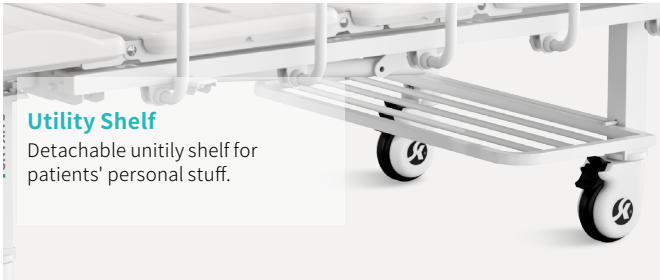

Safe working load

200_{kg}



Back-rest adjustment

70°



Knee-rest adjustment

38°



CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service

WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
 W3400 Corrientes, Argentina
 TEL: +5493794312231
 Av. Senador Vergueiro, 4204 - sala 31
 Rudge Ramos - São Bernardo do Campo - SP, Brazil
 TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
 Mezzanine Floor, Taj Tower,
 upper hill, nairobi, Kenya
 TEL: +254714084806

South East Asia Branch

Saikang Medical Technology Ltd
 Sala Thammasop Subdistrict
 Thawi
 Watthana District
 Bangkok Thailand
 TEL: +66 642741046

UK Branch

Uk Healthward International Limited
 DEPT 111,
 196 High Road Wood Green,
 London, N22 8HH
 United Kingdom
 TEL: +31 (0) 202111106

Russia Branch

LLC 'SAIKANG'
 454004, Russian Federation,
 Chelyabinsk,
 Universitetskaya Naberezhnaya,
 Russia
 TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
 Musaffah
 P.O Box 53200, Abu Dhabi,
 UAE
 Tel: +971 585710652

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
 STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
 📞 +86 18021231901 · 🕒 24/7
 ✉️ export@saikangmedical.com
 🌐 www.saikangmedical.com



Company video



General catalogue